



WE MADE SLEEP THERAPY PRODUCTS EASIER TO USE SO IT'S EASIER FOR YOU TO IMPROVE PATIENTS' LIVES.

When your patients shop for consumer products, they're looking for designs that will **make life easier** – and **fit their lifestyle**. We've found that they want the same qualities in their sleep therapy products. That's why the REMstar M Series offers an incredibly easy-to-understand interface, pressure delivery performance and a sleek design that helps make sleep therapy easier to integrate into a patient's life. And it's another simple example of how Respiration's **big ideas make a real difference** in your business – and your patients' lives.

REMstar® M Series

REMstar® M Series


Sleek, simple CPAP therapy that your patients can trust.

Respironics' REMstar M Series is a full-featured unit that gives patients more of what they want without taking away anything you need. It is smaller, sleeker and easier to use to better fit with a patient's sleep environment and lifestyle. And the new REMstar M Series offers new phone-in compliance reporting capabilities with Encore® Pro for easier patient management.

The Features

- Deluxe ramp capabilities allow a 0–45 minute period with a patient-adjustable ramp start pressure
- Lighted patient buttons control the CPAP unit
- Improved carrying case offers exceptional padding and is big enough to carry the CPAP unit, the mask and the heated humidifier
- Pressure calibration screen allows you to set target pressure to within 0.1 cm H₂O
- 3-position altitude compensation
- Quiet operation

Encore® Pro Compliance Reporting Capabilities

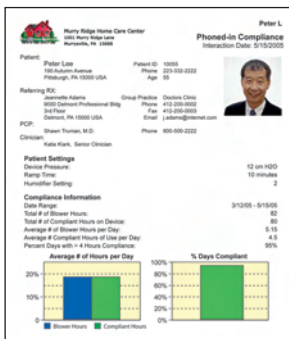

- Data display for phone-in compliance
- Easy report generation featuring:
 - Average hours of use per day
 - % of days compliant

REMstar Heated Humidifier

- Unique design works with all new REMstar M Series sleep therapy devices
- New “hidden” water chamber is easier to use and easier to fill, and offers eight full hours of continuous use
- New easy-turn dial offers precise settings that don't need to be reset if therapy is temporarily suspended



Product Specifications

	REMstar M Series
PRESSURE	4 to 20 cm H ₂ O
RAMP TIME	0 to 45 min. (5-min. increments)
STARTING RAMP PRESSURE	4 to CPAP; patient-adjustable
DIMENSIONS	7.5" L x 5.0" W x 3.125" H 19 cm L x 12.7 cm W x 7.9 cm H
WEIGHT	2.2 lbs
FILTERS	Pollen and optional ultra-fine
DEVICE SET-UP	LCD/keypad
DATA STORAGE CAPACITY (minimum)	Therapy meter
DATA STORAGE CONTENT	Total therapy hours Total blower hours Usage sessions > 4 hrs
COMPLIANCE MONITORING	Breathing detection
ALTITUDE COMPENSATION	Low, medium or high
ELECTRICAL REQUIREMENTS	100 – 240 VAC, 50/60 Hz
HUMIDIFICATION	Integrates with M Series Heated Humidifier and M Series Pass-over Humidifier
DC POWER	Direct connect cord
WARRANTY	2 years
MISC	Fine pressure calibration, lighted LEDs

Ordering Information

REMstar Systems	Part #
REMstar (device only)	DS100
REMstar Core (includes humidifier)	DS100H
CPAP Accessories	Part #
POLLEN FILTERS	1029330
ULTRA-FINE FILTERS	1029331
CARRYING CASE	1022219
TRAVEL PACK	1028122
DC ADAPTER CORD	1001956
SHIELDED DC CORD SYSTEM (incl. shielded DC cord and battery adapter cable)	1001979
REUSABLE PATIENT TUBING (with swivel)	1022333
TUBING SWIVEL	1029532
AC POWER CORD (3 ft / ~1m)	1028078
EXTERNAL POWER SUPPLY	1015642
M SERIES HEATED HUMIDIFIER	1022257
M SERIES PASS-OVER HUMIDIFIER	1029090
M SERIES HUMIDIFIER CHAMBER (replacement)	1029533
ENCORE PRO SOFTWARE	1009583



Customer Service: 800-345-6443 or 724-387-4000

Respironics Europe: +33-1-47-52-30-00

Respironics Asia Pacific: +81-3-5280-9611

www.mseries.respironics.com

CAUTION: U.S. federal law restricts this device to sale by or on the order of a physician.

Respironics, REMstar and Encore are registered trademarks of Respironics, Inc. and its affiliates.

©2006 Respironics, Inc. and its affiliates. All rights reserved.