




New Easy-to-Fit CPAP Interface

The *DreamFit*™ nasal mask from Puritan Bennett expands our selection of CPAP interfaces. This cost-effective mask features simple adjustments that make it easy to fit a wide range of patients. And the open design minimizes contact with the face, so it's comfortable to wear night after night.

[dreamfit™ nasalmask]

- > **The *DreamFit* nasal mask has a number of features to enhance comfort and allow greater freedom of movement.**
 - Airway hose routes over the head to enable sleep in virtually any position.
 - Flexible “bellows” chamber between the mask and interface swivels to accommodate patient movement and help prevent leaks.
 - Soft, pliable mask creates an effective, yet comfortable seal.
 - Quick-release clips make it easy for patients to apply and remove interface.
 - Soft, stretchable Breath-O-Prene® straps keep moisture away from the skin for greater comfort.
 - Vent hole is angled upward to direct air away from bed partner.
 - Comes in three mask sizes (shallow, standard and large) to fit different facial sizes and shapes.

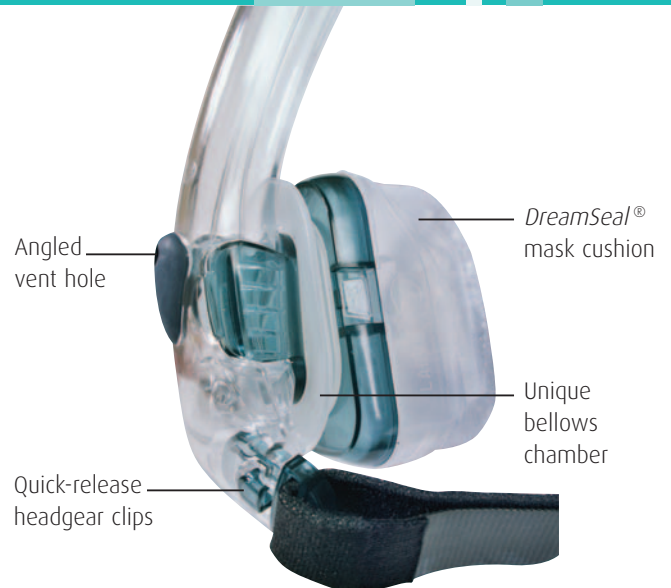

**PURITAN
BENNETT**

[dreamfitnasalmask]

Specifications

Materials of Construction

- > *DreamFit* Mask:
Perimeter seal is medical-grade silicone
- > Straps: Breath-O-Prene
- > All *DreamFit* components are latex-free

Pressure Range

- > 3 – 35 cm H₂O

Hose Attachment

- > Fits 22 mm O.D. connection



DreamFit Nasal Mask

Ordering Information

| Description | Catalog Number |
|--|----------------|
| <i>DreamFit</i> Interface with Large Mask | DF301 |
| <i>DreamFit</i> Interface with Standard Mask | DF302 |
| <i>DreamFit</i> Interface with Shallow Mask | DF303 |

Replacement Parts

| Description | Catalog Number |
|---|----------------|
| <i>DreamSeal</i> Adapter Kit (bellows with clip, <i>DreamSeal</i> adapter) | 10001990 |
| <i>DreamFit</i> Strap Kit (headgear, clips) | 10001994 |
| <i>DreamFit</i> Hose Kit (hose, swivel adapter) | 10002002 |
| <i>DreamFit</i> Silicone Parts Kit (exhaust vent, bellows, forehead cushion supports) | 10002003 |

The *DreamFit* nasal mask is part of our complete family of sleep therapy devices and interfaces, including:

GoodKnight® 420 Series CPAP System

- > Highly responsive to breathing patterns for greater comfort.
- > Amazingly small and lightweight (1.5 lbs).
- > Optional heated humidifier available.



GoodKnight 420 CPAP System

Breeze SleepGear™ CPAP Interface System

- > Over-the-head design allows patients to sleep in virtually any position.
- > Can be used with Puritan Bennett's *DreamSeal* mask or Nasal Pillows, which some patients find to be more comfortable.



Breeze SleepGear Interface With Nasal Pillows

Caution: Federal (U.S.A) law restricts this device to sale by or on the order of a physician.

Specifications subject to change without notice.

*Breath-O-Prene is a registered trademark of AccuMED Technologies, Inc. All other trademarks belong to Tyco Healthcare Group LP or its affiliate. © 2005 Nellcor Puritan Bennett Inc. All rights reserved.

A.ae 2175v2-0905 ST06604

To learn more about our complete line of sleep therapy devices, contact your local Puritan Bennett representative, call 1-800-635-5267 or visit www.puritanbennett.com.

tyco
Healthcare

**Puritan
Bennett**

4280 Hacienda Drive
Pleasanton, CA 94588
Tel 925.463.4000
Toll Free 1.800.635.5267
www.puritanbennett.com

Tyco Healthcare UK LTD.
154 Fareham Road
Gosport, UK PO13 0AS
Tel +44.1329.224000