


# FlexiFit™ 406

PETITE NASAL MASK

FlexiFit™ 407

FlexiFit™ 406

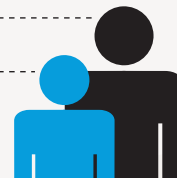

Provided for your reference by

**cpapXchange.com**

## Specification Sheet

A Nasal Delivery Mask




*Specifically designed to fit finer featured faces. Incorporating FlexiFit™ Technology with auto-contouring properties, the FlexiFit™ 406 Nasal Mask comes with a unique Glider™ Strap mechanism for maximum comfort and stability. Other unique features ensure a premium sleep experience.*

**SleepStyle™**  
SLEEP PERFORMANCE FOR  
AN ENERGIZED LIFESTYLE

**Fisher & Paykel**  
HEALTHCARE




# Nasal Delivery - FlexiFit™ 406

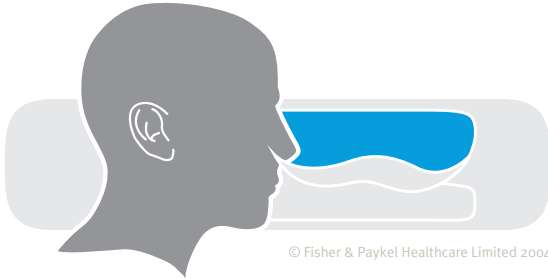

© Fisher & Paykel Healthcare Limited 2004


The FlexiFit™ 406 Petite Nasal Mask offers an over-the-nose mask solution specifically designed to fit finer feature faces. The Petite design is small, light and easy to use. The Stretchgear™ Headgear provides even greater comfort and a maintenance-free bias diffuser ensures quieter operation. The FlexiFit™ 406 Petite Nasal Mask comes as an Out-of-Box Solution fully assembled and ready for use.

## UNIQUE MASK FEATURES FOR GREATER FREEDOM, COMFORT AND COMPLIANCE

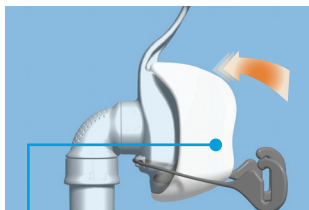

### FlexiFit™ Technology Auto-contouring fit

*With Foam Cushion*  
FlexiFit™ Technology automatically adjusts to a wide range of nasal and facial contours with simple one-step ease of fitting and no complicated adjustments.

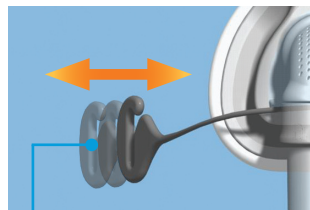

### Glider™ Strap Anti-leak system

The sliding Glider™ Strap provides freedom of movement while maintaining a seal. Mask force remains equal when the head rotates from side to side, minimizing the occurrence of leaks.




### Stretchgear™ Headgear Deluxe mask suspension system

This is a special elastic headgear segment that maintains an equal mask force when the head tilts up or down.




### Premium Frosted Silicone Seal Optimizing seal technology

This seal technology has a soft, smooth feel to slide onto the face, and its frosted texture optimizes the mask seal.

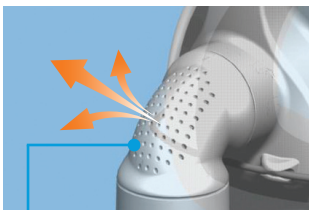


### Efficient Bias Diffuser Quiet usage

*Maintenance Free*  
This advanced air diffuser system prevents air from blowing onto the sleeping partner and decreases noise that could potentially disrupt sleep.

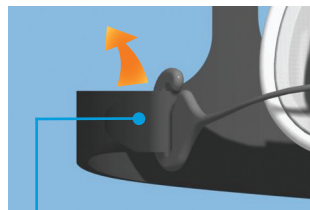

### Glider™ Clip One-step removal

Simple one-step removal of headgear. No ongoing headgear adjustments required.

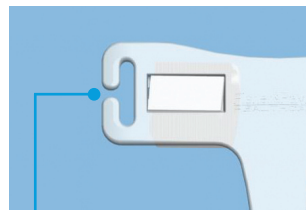

### Hidden Strap Release One-step removal

Easy one-step removal of headgear. Enables reassembly without having to make fitting adjustments.

## CONTENTS

- Complete mask with headgear assembled and ready for use
- Single one size Silicone Seal designed to fit most patients
- User friendly Fitting Guide

REORDER CODE: HC406

FOR MORE INFORMATION, PLEASE CONTACT YOUR LOCAL FISHER AND PAYKEL HEALTHCARE REPRESENTATIVE:

**International** PO Box 14 348, Panmure, Auckland 1134, New Zealand Tel: +64 9 574 0100 Fax: +64 9 574 0158 Email: info@fphcare.com Web: www.fphcare.com

**USA** Tel: 1800 446 3908 +1 949 470 3900 Fax: +1 949 470 3933 **Australia** Tel: +61 3 9879 5022 Fax: +61 3 9879 5232 **Austria** Tel: 0800 29 31 23 Fax: 0800 29 31 22 **Benelux** Tel: +31 78 644 0924 Fax: +31 78 644 0914 **China** Tel: +86 20 8363 4286 Fax: +86 20 8363 3595 **France** Tel: +33 1 6446 5201 Fax: +33 1 6446 5221 **Germany** Tel: +49 7182 93777 0 Fax: +49 7182 93777 99 **India** Tel: +91 80 5123 6041 Fax: +91 80 5123 6044 **Irish Republic** Tel: 1800 409 011 **Italy** Tel: +39 06 7839 2939 Fax: +39 06 7814 7709 **Japan** Tel: +81 3 3661 7205 Fax: +81 3 3661 7206 **Spain** Tel: +34 902 013 346 Fax: +34 902 013 379 **Switzerland** Tel: 0800 83 47 63 Fax: 0800 83 47 54 **UK (EU Authorised Representative)** Tel: +44 1628 626 136 Fax: +44 1628 626 146 **Nthn Ireland** Tel: 0800 132 189