

**Women everywhere will thank
you in the morning.**



Mojo™...the most comfortable SleepMask™ in the world.

Dual beam element transfers strap forces into the mask for an exceptionally secure fit

Forehead pad stabilizes mask with soft gel cushion


Proprietary tensioning screw provides extremely fine adjustment for perfect sealing

Gel cushion yields to soft face tissue

Magnetic clips provide easy capture and secure fastening

Non-rebreathing valve provides an integral safety feature in case of PAP device failure

Pliable ring and flexible shell for instant customization

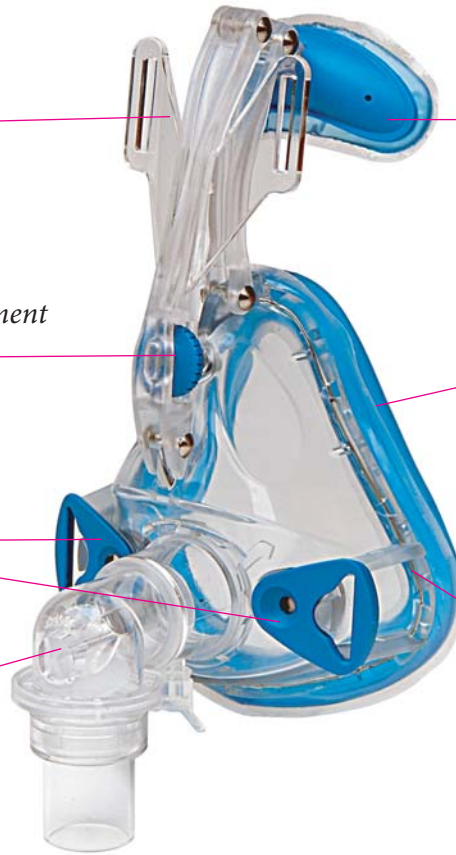

More comfortable, more sellable.

The superb comfort of the Mojo™ SleepMask™ is a result of combining physics and engineering with state-of-the-art materials.

Mojo's contoured shell follows the curves of the human face, and easily bends inward and outward for a comfortable fit. Air leakage is eliminated by adjusting the headgear, as the dual beam element transfers the force of the straps to the center of the mask. Our proprietary tension wheel fine-tunes the fit.

SleepNet's unique magnetic clips are virtually foolproof. They easily click into position securing the mask and keeping the straps in place. The clips also allow for a quick release of the mask. The SleepNet Full Face Mask, designed using strong, durable and medically approved materials, will easily withstand the rigors of lab or home use. Best of all, the outstanding features and comfort of the Mojo SleepMask will meet the needs and demands of both patients and technicians.