

# Oracle™ 452

ORAL MASK

## Specification Sheet


---

### An Oral Delivery Mask



*The unique Oracle™ 452 delivers pressure through the mouth rather than the nose. Clinically proven to function as effectively as a nasal mask\*, the Oracle™ 452 is a perfect solution for mouth breathers and patients with chronic nasal obstruction.*

# Oral Delivery - Oracle™ 452

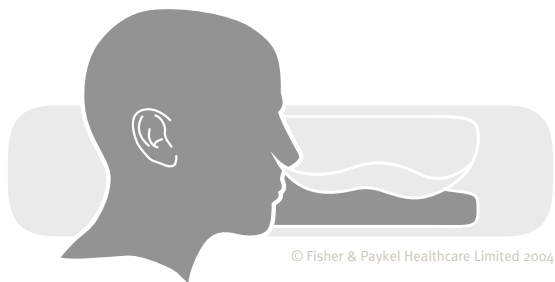

*The Oracle™ 452 Oral Mask is easy to use and easy to seal. With nothing covering the field of vision, the Oracle™ 452 provides freedom of movement equivalent only to that of the Infinity™ 481 Direct Nasal Mask. The Oracle™ 452 Oral Mask comes as an Out-of-Box Solution fully assembled and ready for use.*

## UNIQUE MASK FEATURES FOR GREATER FREEDOM, COMFORT AND COMPLIANCE




### Single Strap Headgear Minimal Intrusion

Minimal, lightweight and comfortable, single BREATH-O-PRENE™ strap allows clear vision and movement. Patient has freedom to choose Flexitube position with 360° of movement.




### Rotating Adjustment Dial Easy Adjustment

Offers continuous adjustment of the distance between the SnapFlap™ Cover and SoftSeal Mouthpiece for a personalized fit. Easy to use - just rotate the adjustment dial. This can even be adjusted while mask is fitted.

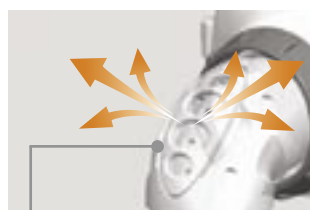

### Efficient Bias Diffuser Quiet usage

**Material Filter**  
This advanced air diffuser system prevents air from blowing onto the sleeping partner and decreases noise that could potentially disrupt sleep.

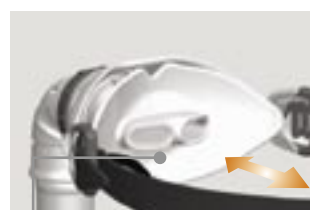

### SoftSeal Mouthpiece/ SnapFlap™ Cover Interaction

Holds the lips snugly between the inner SoftSeal Mouthpiece and outer SnapFlap™ Cover providing a leak-free combination. These technologies work together to provide a simple, superior seal and maximum comfort.



### Easy-Release Strap One-step removal

Simply pull the Easy-Release Cord for one-step removal of headgear. No ongoing headgear adjustments required.

## CONTENTS

- Assembled mask complete with Headgear and Nasal Plugs
- Alternate size Silicone Softseal Mouthpiece included (S,L)
- User friendly Fitting Guide

REORDER CODE: HC452

## REFERENCES

\* Smith PL, O'Donnell CP, Allan L et al. A physiologic comparison of nasal and oral positive airway pressure. *Chest* 2003; 123 (3):689-94. \* Anderson FE, Kingshott RN, Taylor DR et al. A randomized crossover efficacy trial of oral CPAP (Oracle) compared with nasal CPAP in the management of obstructive sleep apnea. *Sleep* 2003; 26 (6):721-6. \* Khanna R, Kline LR. A prospective 8 week trial of nasal interfaces vs. a novel oral interface (Oracle) for treatment of obstructive sleep apnea hypopnea syndrome. *Sleep Med* 2003; 4 (4):333-8. \* Beecroft J, Zanon S, Lukic D et al. Oral continuous positive airway pressure for sleep apnea: effectiveness, patient preference, and adherence. *Chest* 2003; 124 (6):2200-8.

FOR MORE INFORMATION, PLEASE CONTACT YOUR LOCAL FISHER AND PAYKEL HEALTHCARE REPRESENTATIVE:

**International** P O Box 14 348, Panmure, Auckland 1134, New Zealand Tel: +64 9 574 0100 Fax: +64 9 574 0158 Email: info@fphcare.com Web: www.fphcare.com

**USA** Tel: 1800 446 3908 +1 949 470 3900 Fax: +1 949 470 3933 **Australia** Tel: +61 3 9879 5022 Fax: +61 3 9879 5232 **Austria** Tel: 0800 29 31 23 Fax: 0800 29 31 22 **Benelux** Tel: +31 78 644 0924 Fax: +31 78 644 0914 **China** Tel: +86 20 8363 4286 Fax: +86 20 8363 3595 **France** Tel: +33 1 6446 5201 Fax: +33 1 6446 5221 **Germany** Tel: +49 7182 93777 0 Fax: +49 7182 93777 99 **India** Tel: +91 80 5123 6041 Fax: +91 80 5123 6044 **Irish Republic** Tel: 1800 409 011 **Italy** Tel: +39 06 7839 2939 Fax: +39 06 7814 7709 **Japan** Tel: +81 3 3661 7205 Fax: +81 3 3661 7206 **Spain** Tel: +34 902 013 346 Fax: +34 902 013 379 **Switzerland** Tel: 0800 83 47 63 Fax: 0800 83 47 54 **UK (EU Authorised Representative)** Tel: +44 1628 626 136 Fax: +44 1628 626 146 **Nthn Ireland** Tel: 0800 132 189