


F&P Simplus™

Simply Fits + Performs

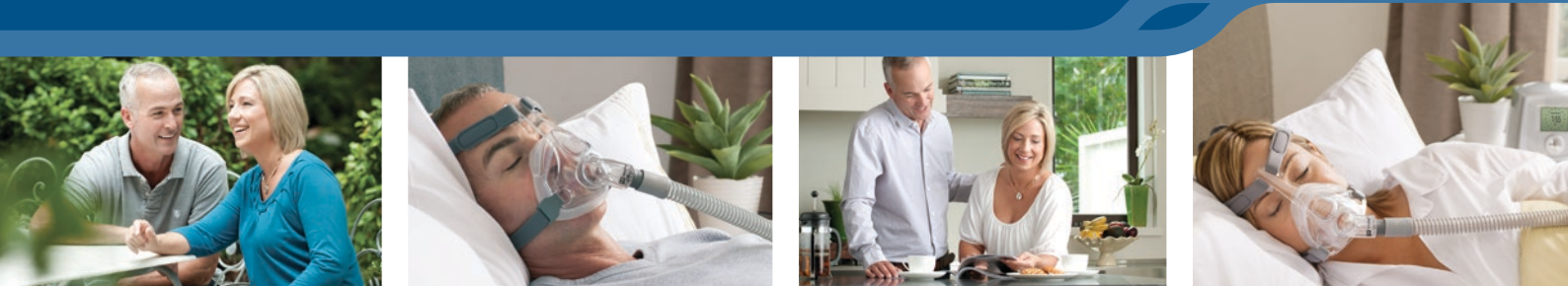

F&P SIMPLUS FULL FACE MASK

When engineering the F&P Simplus our designers set out to create a mask that revolutionized Full Face comfort, seal and easy use.

The F&P Simplus is the result of intensive research and development, many prototypes and thousands of design hours. We consulted patients, equipment providers and clinicians from around the world to develop a mask that **simply fits + performs.**

F&P Simplus Features

The F&P Simplus incorporates three key components, the **RollFit™ Seal**, **ErgoForm™ Headgear** and **Easy Frame**, all designed to work in harmony. In combination, these components offer the comfort, seal and easy use that Fisher & Paykel Healthcare masks are known for.

ErgoForm Headgear
Easy-Clip Hooks +
Adjustable Crown Strap +
Stretch and Non-stretch Panels +
Embedded VELCRO® Tabs +

RollFit Seal
Integrated Stability Panels +
Advanced Air Diffuser +

Easy Frame
Easy-Clip Frame Attachment +
Ball-and-Socket Elbow +



1. RollFit Seal

As the name suggests, the one-piece seal “rolls” back and forth on the bridge of the nose to adjust automatically. In doing so, the **RollFit Technology** reduces pressure on the bridge of the nose.

The **RollFit Seal** incorporates **Integrated Stability Panels** to enhance the sealing performance and the **Advanced Air Diffuser** is designed for minimal noise and draft.

2. ErgoForm Headgear

This breathable, ergonomically designed headgear self-locates high on the rear of the head, allowing for unrestricted head movement (sideways, up and down).

Incorporating **Stretch and Non-stretch Panels**, the **ErgoForm Headgear** is designed to hold its shape and allow stretch where it is required, for optimal usability and performance.

Embedded VELCRO®* Tabs with **Easy-Clip Hooks** enable fast fitting and removal.

The **Adjustable Crown Strap** enhances mask stability.

3. Easy Frame

This low-profile frame is stable, durable and small, and ensures a clear line of sight.

The one frame fits all three seal sizes and has an **Easy-Clip Frame Attachment** to aid assembly after cleaning.

The **Ball-and-Socket Elbow** rotates for freedom of movement while reducing drag from the CPAP tube.

Mask Reorder Codes

| | |
|---------------|---------------------------|
| 400475 | F&P Simplus Mask – Small |
| 400476 | F&P Simplus Mask – Medium |
| 400477 | F&P Simplus Mask – Large |

Spare Parts Reorder Codes

| | | | |
|-----------------|---------------------------------------|-----------------|-------------------------------------|
| 400HC579 | F&P Simplus Seal – Small | 400HC582 | F&P Simplus Headgear – Small |
| 400HC580 | F&P Simplus Seal – Medium | 400HC583 | F&P Simplus Headgear – Medium/Large |
| 400HC581 | F&P Simplus Seal – Large | 400HC584 | F&P Simplus Mask No H/G – Small** |
| 400HC570 | F&P Simplus Swivel | 400HC585 | F&P Simplus Mask No H/G – Medium** |
| 400HC587 | F&P Simplus Headgear – Clips & Buckle | 400HC586 | F&P Simplus Mask No H/G – Large** |

* VELCRO® is a registered trademark of Velcro Industries B.V.
** Not available in all countries