



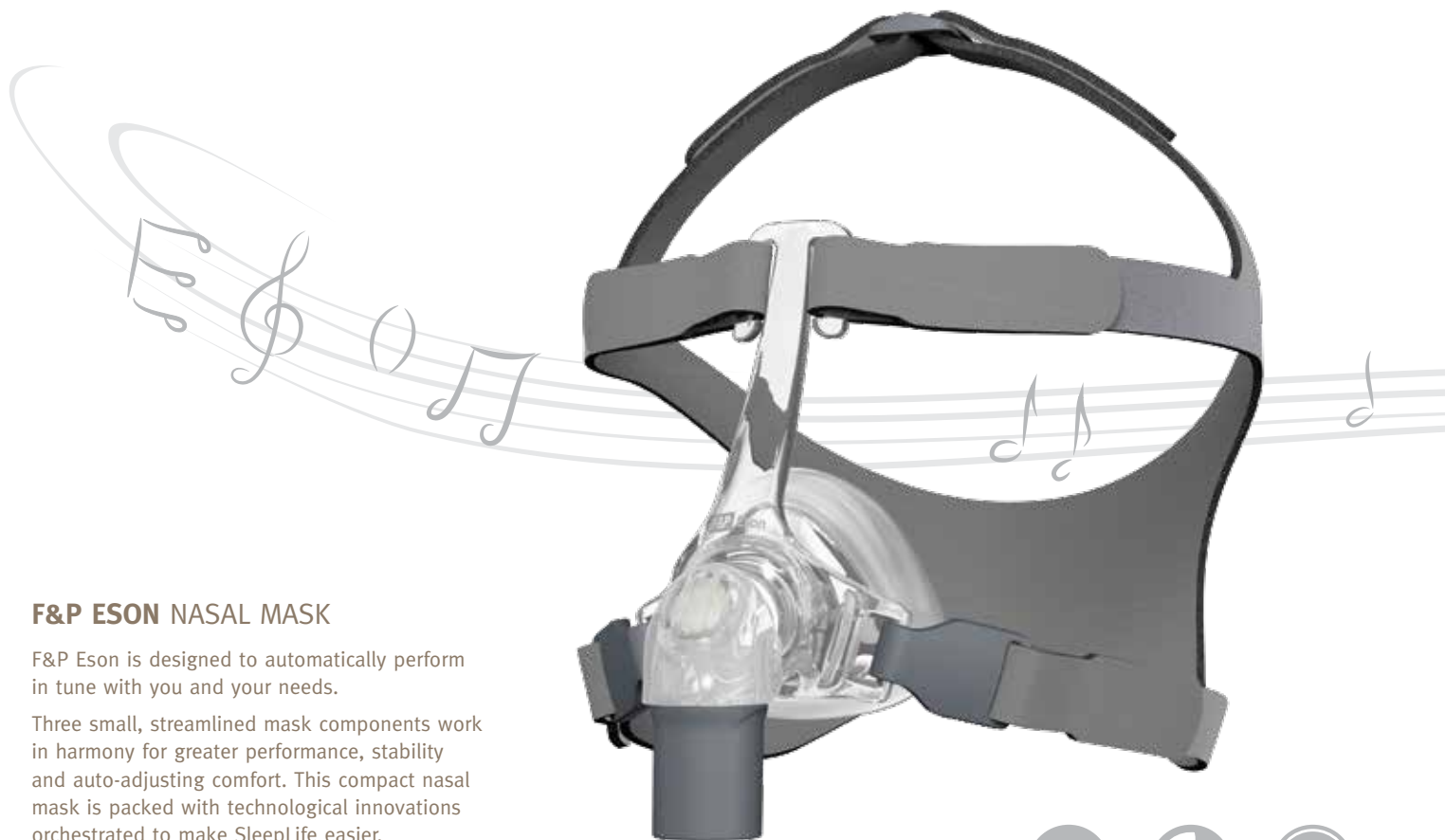

## Performs Automatically - in Tune with You

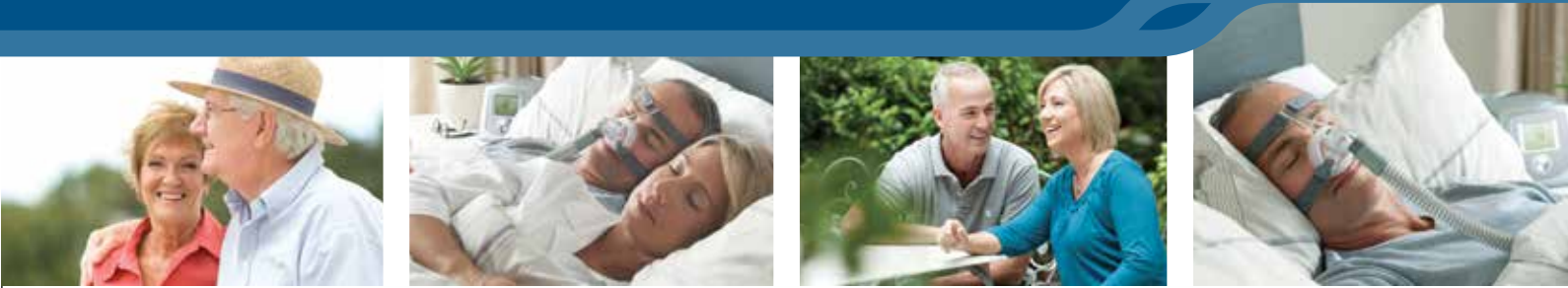

### F&P ESON NASAL MASK

F&P Eson is designed to automatically perform in tune with you and your needs.

Three small, streamlined mask components work in harmony for greater performance, stability and auto-adjusting comfort. This compact nasal mask is packed with technological innovations orchestrated to make SleepLife easier.



Comfort



Seal



Easy Use

F&P Eson is very 'easy on you' from bag to bed. It has been designed to perform automatically in tune with every aspect of SleepLife:\* F&P Eson maximizes no-fuss comfort and simplicity by minimizing manual adjustments, parts and noise.

## EASY ON PROVIDERS AND SLEEPLABS

## EASY ON PATIENTS AND THEIR PARTNERS

### Easy Fit



Intuitive fitting innovations make it easy and fast to fit your patients.

Easy-Clip hooks allow for fast fitting and removal.

With only three simple parts, it is easy to clean and reassemble.


Design innovations maximize comfort and optimize a stable, effective seal all night.

It is super-quiet with minimal draft allowing an undisturbed sleep for you and your partner.



\* We sleep for a third of our lives. During sleep, we grow, we repair and memories are made; our bodies and minds are rejuvenated. Sleep feeds our quality of life.

F&P Eson's three simple components, the **RollFit™ Seal**, **ErgoFit Headgear** and **Easy Frame**, work in harmony to deliver the comfort, seal and easy use that Fisher & Paykel Healthcare masks are known for.

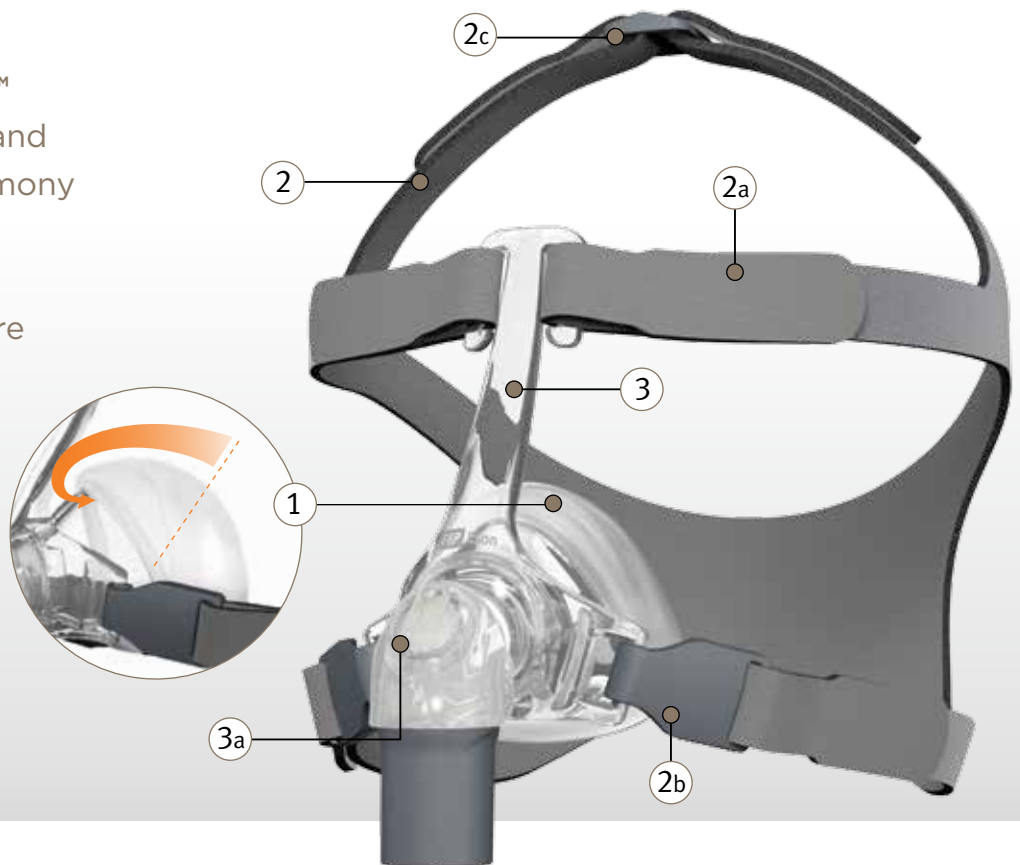

### 1. RollFit Seal

*Easy auto-adjusting fit that optimizes an effective, comfortable seal*

As the name suggests, the one-piece seal 'rolls' back and forth on the bridge of the nose to adjust automatically. In doing so, the **RollFit Technology** minimizes pressure on the bridge of the nose.

The **RollFit Seal** incorporates **Integrated Stability Panels** for optimal sealing performance and stability all night long.

### 2. ErgoFit Headgear

*Easy to use, fit, adjust and clean*

This breathable, ergonomically designed headgear has been designed to self-locate high on the rear of the head allowing for maximum head movement (sideways, up and down) without mask dislodgement.

The **ErgoFit Headgear** also incorporates **Embedded VELCRO®\*** **Tabs** (2a) with **Easy-Clip hooks** (2b) for fast and easy fitting and removal.

The adjustable **Crown Strap** (2c) enhances mask stability.

\* VELCRO® is a registered trademark of Velcro Industries B.V.

### 3. Easy Frame

*Easy to wear, see and sleep*

This low-profile frame is stable, durable and small, and ensures a clear line of sight.

The one frame fits all three seal sizes and has an **'Easy-Clip'** frame attachment for effortless assembly after cleaning.

The **Easy Frame** comes with the **Q Cover** (3a) for super-quiet use and minimal partner draft.

The Ball and Socket Elbow rotates for freedom of movement while minimizing drag from the CPAP tube.

F&P Eson performs in tune with what patients and their healthcare providers want: **Automatic Performance** that makes adherence **Easy on Everyone**.

For more information, please contact  
your local Fisher & Paykel Healthcare representative

**Manufacturer** 

Fisher & Paykel Healthcare Ltd  
15 Maurice Paykel Place  
East Tamaki, Auckland 2013

PO Box 14 348 Panmure  
Auckland 1741  
New Zealand

Tel: +64 9 574 0100  
Fax: +64 9 574 0158  
Email: [info@fphcare.co.nz](mailto:info@fphcare.co.nz)  
Web: [www.fphcare.com](http://www.fphcare.com)

**Australia**

Fisher & Paykel Healthcare Pty Limited  
36-40 New Street, PO Box 167  
Ringwood, Melbourne  
Victoria 3134, Australia

Tel: +61 3 9879 5022  
Fax: +61 3 9879 5232

**Austria**

Tel: 0800 29 31 23  
Fax: 0800 29 31 22

**Benelux**

Tel: +31 40 216 3555  
Fax: +31 40 216 3554

**China**

Tel: +86 20 3205 3486  
Fax: +86 20 3205 2132

**France**

Tel: +33 1 6446 5201  
Fax: +33 1 6446 5221

**Germany**

Tel: +49 7181 98599 0  
Fax: +49 7181 98599 66

**India**

Tel: +91 80 4284 4000  
Fax: +91 80 4123 6044

**Irish Republic**

Tel: 1800 409 011

**Italy**

Tel: +39 06 7839 2939  
Fax: +39 06 7814 7709

**Japan**

Tel: +81 3 3661 7205  
Fax: +81 3 3661 7206

**Northern Ireland**

Tel: 0800 132 189

**Spain**

Tel: +34 902 013 346  
Fax: +34 902 013 379

**Sweden**

Tel: +46 8 564 76 680  
Fax: +46 8 36 63 10

**Switzerland**

Tel: 0800 83 47 63  
Fax: 0800 83 47 54

**Taiwan**

Tel: +886 2 8751 1739  
Fax: +886 2 8751 5625

**Turkey**

Tel: +90 312 354 34 12  
Fax: +90 312 354 31 01

**UK**



Fisher & Paykel Healthcare Ltd  
Unit 16, Cordwallis Park  
Clivemont Road,  
Maidenhead  
Berkshire SL6 7BU, UK

Tel: +44 1628 626 136  
Fax: +44 1628 626 146

**USA/CANADA**

Tel: +1 800 446 3908  
or +1 949 453 4000  
Fax: +1 949 453 4001

フィッシャー & パイケルヘルスケア株式会社  
〒103-0025 東京都中央区日本橋茅場町3-8-10 リベラ茅場町5F  
FAX: 03-3661-7206 ㊦注文書専用 FAX: 0800-100-7207  
TEL: 03-3661-7205 (代表)

Mask Reorder Codes		Spare Parts Reorder Codes			
<b>400449</b>	F&P Eson Mask Small	<b>400HC567</b>	F&P Eson Headgear - Small	<b>400HC228</b>	F&P Eson Diffuser (10pk) and Cover
<b>400450</b>	F&P Eson Mask Medium	<b>400HC568</b>	F&P Eson Headgear - Medium/Large	<b>400HC570</b>	F&P Eson Swivel
<b>400451</b>	F&P Eson Mask Large	<b>400HC569</b>	F&P Eson Headgear - Clips & Buckle	<b>400HC571</b>	F&P Eson Mask No H/G - Small
		<b>400HC126</b>	F&P Eson Seal - Small	<b>400HC572</b>	F&P Eson Mask No H/G - Medium
		<b>400HC127</b>	F&P Eson Seal - Medium	<b>400HC573</b>	F&P Eson Mask No H/G - Large
		<b>400HC128</b>	F&P Eson Seal - Large		