

FlexiFit™ 432

FULL FACE MASK

Provided for your reference by
cpapXchange.com


Specification Sheet

A Full Face Delivery Mask




Our second generation Full Face Mask incorporates a soft Foam Cushion that works together with the highly contoured Silicone Seal to provide improved sealing technology and a more comfortable fit.

Full Face Delivery - FlexiFit™ 432

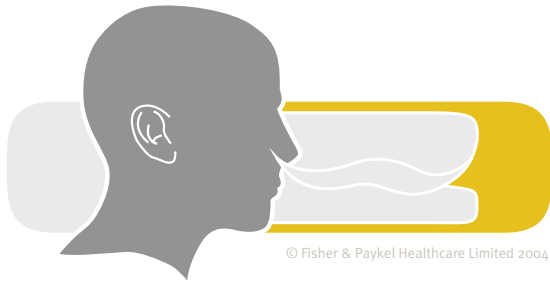


The FlexiFit™ 432 Full Face Mask offers an improved seal and more comfortable fit with the inclusion of a soft Foam Cushion that interacts with the Silicone Seal. This alternative to nasal or oral delivery comes fully assembled and ready for use.

UNIQU MASK FEATURES FOR GREATER FREEDOM, COMFORT AND COMPLIANCE

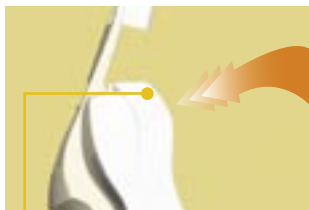


FlexiFit™ Technology Auto-contouring fit

With Foam Cushion FlexiFit™ Technology automatically adjusts to a wide range of nasal and facial contours with simple one-step ease of fitting and no complicated adjustments. The Foam Cushion improves the sealing technology and comfort.

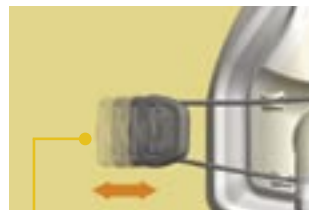

Glider™ Strap Anti-leak system

The sliding Glider™ Strap provides freedom of movement while maintaining a seal. Mask force remains equal when the head rotates from side to side, minimizing the occurrence of leaks.

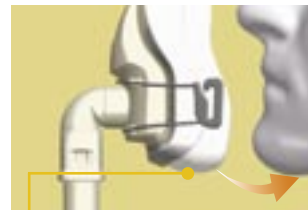

Under-Chin Design Extra stability

An Under-Chin fit provides stability by providing a constant point of reference and maintains the seal even with movement of the jaw during the night.



Efficient Bias Diffuser Quiet usage

Maintenance Free
This advanced air diffuser system prevents air from blowing onto the sleeping partner and decreases noise that could potentially disrupt sleep.




Premium Frosted Silicone Seal Optimizing seal technology

This seal technology has a soft, smooth feel to slide onto the face, and its frosted texture optimizes the mask seal.



Crown Strap Extra stability

Complements the Under-Chin Design to maximize stability and eliminate the need for a chin strap.

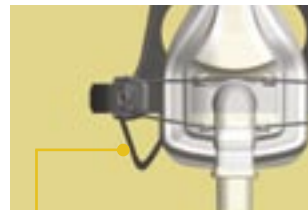

Easy-Release Strap One-step removal

Simply pull the Easy-Release Cord for one-step removal of headgear. No ongoing headgear adjustments required.

CONTENTS

- Complete mask with headgear assembled and ready for use
- The FlexiFit™ 432 Full Face comes in three separate sizes, Small, Medium and Large
- User friendly Fitting Guide

REORDER CODE: HC432

FOR MORE INFORMATION, PLEASE CONTACT YOUR LOCAL FISHER AND PAYKEL HEALTHCARE REPRESENTATIVE:

International P O Box 14 348, Panmure, Auckland 1134, New Zealand Tel: +64 9 574 0100 Fax: +64 9 574 0158 Email: info@fphcare.com Web: www.fphcare.com

USA Tel: 1800 446 3908 +1 949 470 3900 Fax: +1 949 470 3933 **Australia** Tel: +61 3 9879 5022 Fax: +61 3 9879 5232 **Austria** Tel: 0800 29 31 23 Fax: 0800 29 31 22
China Tel: +86 20 3205 3486 Fax: +86 20 3205 3495 **France** Tel: +33 1 6446 5201 Fax: +33 1 6446 5221 **Germany** Tel: +49 7182 93777 0 Fax: +49 7182 93777 99
India Tel: +91 80 4123 6041 Fax: +91 80 4123 6044 **Irish Republic** Tel: 1800 409 011 **Italy** Tel: +39 06 7839 2939 Fax: +39 06 7814 7709 **Japan** Tel: +81 3 3661 7205
UK (EU Authorised Representative) Tel: +44 1628 626 136 Fax: +44 1628 626 146 **Nthn Ireland** Tel: 0800 132 189
Netherlands Tel: +31 78 644 0924 Fax: +31 78 644 0914 **Spain** Tel: +34 902 013 346 Fax: +34 902 013 379 **Switzerland** Tel: 0800 83 47 63
USA Fax: +81 3 3661 7206