

Provided for your reference by
cpapXchange.com



Right on the nose

Respironics Simplicity Nasal Mask bridges the gap between performance and comfort.

Sometimes achieving effective therapy can compromise comfort. So we developed Simplicity to eliminate pressure points, while still delivering high performance.

Simplicity is a small, lightweight nasal mask worn just on the tip of the nose. By having no contact with the bridge of the nose, it helps avoid common pressure points. As a result, both men and women wearers can enjoy a positive seal without irritation. Patients also will enjoy an unobstructed view – in fact, they can even wear their eyeglasses. Simplicity is one more way Philips Respironics is providing sensible solutions that help make it easier for you to fit 100% of your patients, 100% of the time.

PHILIPS
RESPIRONICS


The exhalation port results in quiet air flow, and its innovative built-in design eliminates the need for additional valves or ports.

360° swivels at the top and bottom allow tubing to twist freely as the patient moves.

The Deluxe Softcap headgear is a lightweight, streamlined mesh cap with adjustable straps that allow for optimal fit.



The Simplicity headgear is specially designed for a secure, comfortable fit. Made of soft, flexible material, it features a unique five-strap design that evenly distributes tension, adds stability, and helps to prevent mask leaks. Velcro tabs allow easy adjustment for a custom fit.

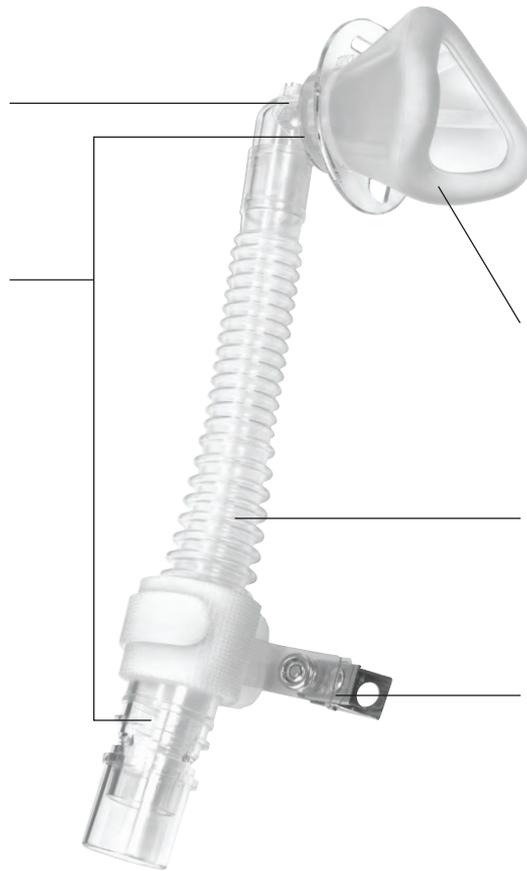

The mask cushion is made of soft, pliable silicone that gently contours around the tip of the nose to provide a snug seal.

Flexible, lightweight tubing is permanently attached to the mask and connects easily to standard circuit tubing.

The Stability Clip on the tubing may be attached to sleepwear or the headgear for added support. A Velcro tab secures placement of the clip on the tubing.

Item	Part number	HCPCS
Simplicity mask		
Mask with Simplicity headgear, small	1002757	A7034, A7035
Mask with Simplicity headgear, medium	1002759	A7034, A7035
Mask with Deluxe Softcap, small	1018138	A7034, A7035
Mask with Deluxe Softcap, medium	1018158	A7034, A7035
Parts and accessories		
Simplicity headgear	1002764	A7035
Deluxe Softcap	1006416	A7035
Replacement cushion, small	1003658	A7032
Replacement cushion, medium	1003659	A7032
Tubing assembly with Stability Clip	1004209	N/A

Respironics, Simplicity, and Softcap are trademarks of Respironics, Inc. and its affiliates, all rights reserved. Velcro is a trademark of Velcro Industries.



© 2009 Koninklijke Philips Electronics N.V. All rights are reserved.

Philips Healthcare reserves the right to make changes in specifications and/or to discontinue any product at any time without notice or obligation and will not be liable for any consequences resulting from the use of this publication.

CAUTION: US federal law restricts these devices to sale by or on the order of a physician.

Hoech JJ 10/19/09 MCI 4102597 PN 1025656

Philips Healthcare is part of Royal Philips Electronics

Respironics Asia Pacific
+65 6298 1088
Respironics Australia
+61 (2) 9666 4444
Respironics Europe, Middle East, Africa
+33 1 47 52 30 00
Respironics United Kingdom
+44 800 1300 845

Philips Respironics
1010 Murry Ridge Lane
Murrysville, PA 15668

Customer Service
+1 724 387 4000
800 345 6443 (toll free, US only)
www.philips.com/respironics